


Optelec TwinView

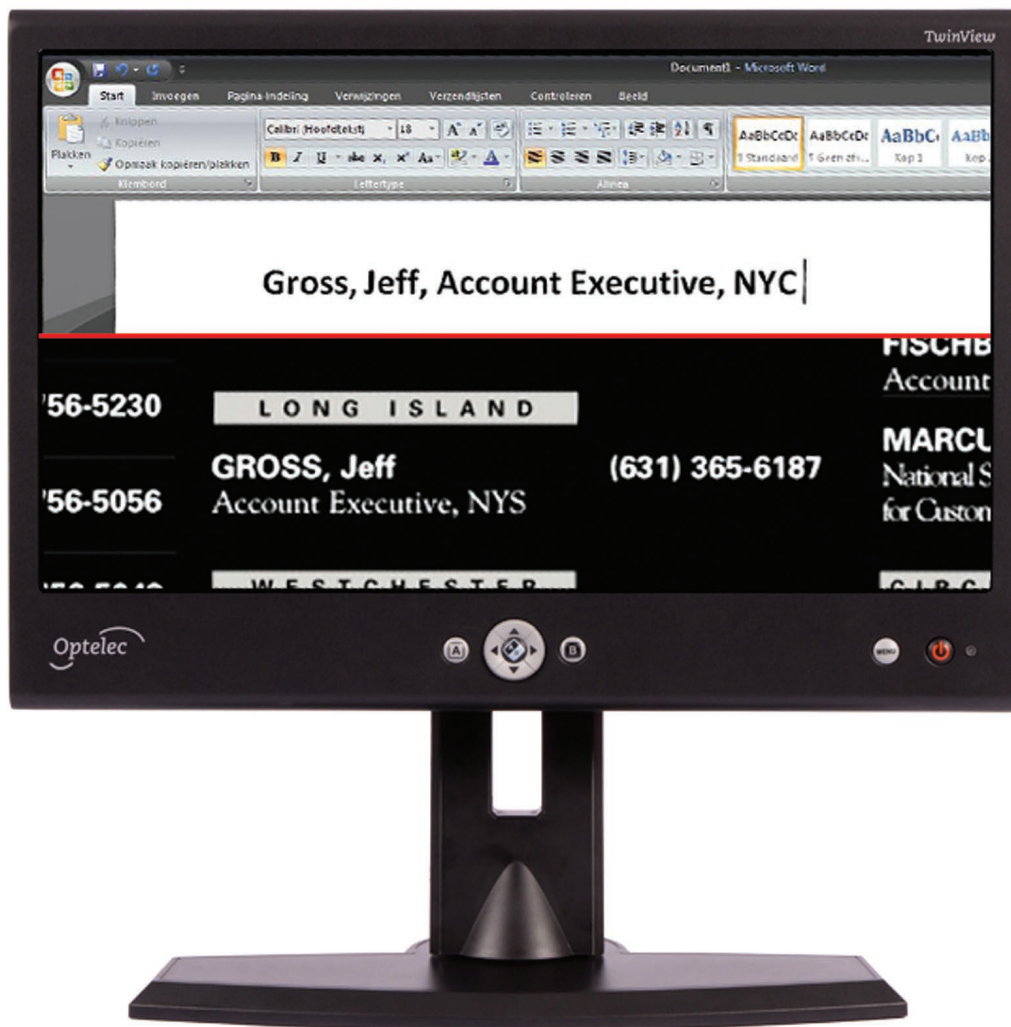

The first Low Vision monitor featuring intelligent split-screen functionality

Imagine using a single low vision workstation, bringing together the functionality of the PC in combination with the facility to reference magnified textbooks and paper documents.

Choose the Optelec TwinView

The Optelec TwinView low vision monitor eliminates the need for multiple screens, harnessing the power of the latest powerful video magnifiers.

Connect the PC and video magnifier to the Optelec TwinView monitor and view images from either device in full screen mode, toggling between both at the press of a button.

Alternatively, using the **UNIQUE** 'True Split-Screen' function, combine and view both images from the PC and video magnifier by choosing one of four configurations and never lose sight of the mouse pointer.

Viewing a PC image and a printed document magnified by a video magnifier on a single screen enables those with

low vision to read and interpret data from printed documents when writing a report, or completing online tax programs whilst viewing paper receipts.

Amazingly, the video magnifier no longer needs to have a special split-screen camera, which is often a more expensive option. The Optelec TwinView low vision monitor provides the split-screen capability, making all the adjustments automatically, and when connected to a high definition video magnifier, ensures the best possible image is displayed at a superior resolution in either full or split-screen.

Four split-screen configurations to choose from:



TwinView

Enjoy true split-screen in glorious HD and never

Add an Optelec Video Magnifier

TwinView low vision screens are available as part of a complete package with any of the Optelec ClearView+ range of video magnifiers and include a Basic Feature Pack offering a range of semi-colours and a position locator:

- TwinView 24-inch Monochrome
- TwinView 24-inch Colour
- TwinView 24-inch HD Colour

Alternatively, connect the TwinView to the Optelec FarView, Sentry or ClearNote video magnifiers, or the ClearReader+ audio text-reader.

TwinView mounted on the ergonomic desktop stand, positioned between the PC and video magnifier



lose sight of the mouse pointer

Optelec
Life is worth enjoying

Choose the Optelec TwinView


Screen

24-inch TFT
High Definition
widescreen
178 degree wide
viewing angle




Connectivity

View full screen (PC or video magnifier)
True split-screen (PC and video magnifier)
Display both images simultaneously using
either horizontal or vertical true split-screen

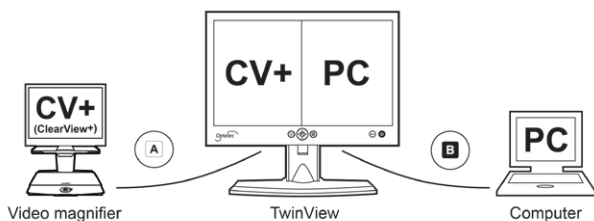

TwinView software plus USB, 2x VGA, 2x DVI and mains power cables included.

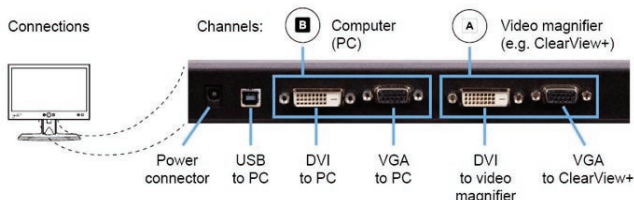

Weight

- 6.7 kg / 14.8 lbs (including desktop stand)
- 4 kg / 8.8 lbs

Designed for low vision

Recognisable buttons
Large font menu
Responsive split-screen
at 60 frames per second



Mounting options

Free-standing, height-adjustable
desktop stand (included)

Position on the
Optelec ClearView+
Video Magnifier

Use a proprietary
monitor arm
(available separately)



Compatibility

- SuperNova (v.12.02 minimum)
- ZoomText (v.9.4 minimum)

Intelligent, fully automatic split-screen
functionality. No further adjustments are
required within the PC magnification
software thanks to the TwinView software.

High specification

- High contrast ratio: 3000:1
- High brightness: 300 cd/m²
- HD resolution: 1920 x 1080 pixels

Freephone 080 8090 8090

to request more information and/or a demonstration:



Optelec Limited, 2 Millfield House, Woodshots Meadow,
Croxley Green Business Park, Watford WD18 8YX.

Phone: 01923 23 13 13 Fax: 01923 23 13 85

info@optelec.co.uk Buy online: www.optelec.co.uk

*“Improving the quality of life of visually impaired and dyslexic people –
reaching out with simple and effective solutions”*